



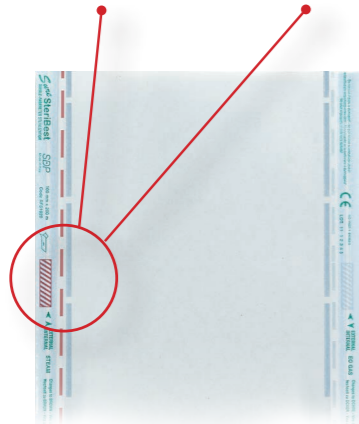

STERILIZATION INDICATORS

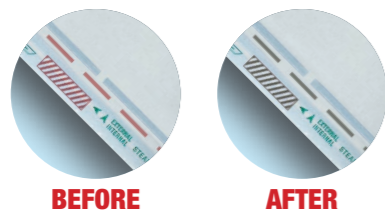

Sure **SteriBest**

FLAT REELS WITH INTERNAL AND EXTERNAL INDICATORS

The special feature that distinguishes these reels from the standard flat reels, is the fact of having in addition to the usual indicator outside, also an indicator within the welding.

The indicator inside is protected by a further welding so that is avoided the contact between the material inserted in the roll and the ink as it is required by the regulations. Benefits of the double parameter of sterilization:

- **Increased safety** during the sterilization process: check done by two indicators, external and internal.
- **Complete guarantee** on the sterilized product: the indicators are in direct contact with the gas that sterilizes the product inserted
- **Saving of money and time:** it is not necessary to insert an indicator strip of sterilization inside the reel/pouch as further check
- **Total and easier check:** regarding the reels, it is possible to check the internal sterilization on all the length of the piece cut
- **Indicators protection:** the internal indicators are protected; they can't enter in contact with the product inserted



A special film with characteristics of resistance and good peelability, sealed to a medical paper with high anti-bacterial protection and good quality to sterilization, are the main components of our flat reels. Raw materials are selected by us, and all the working phases, till to have the finished products, are made under our own direct control in our factory.

FEATURES:

- Medical Paper 60 g/m²
- Laminated Film PET/PP peel blue color
- Indicators STEAM and EO (water based, not toxic)
- Conformity with international standards
- Multi-line weld lines are very durable and equally spaced.
- No release of fibers of paper or lacerations of the plastic material
- Opening direction is marked
- Gusseted Reels allow the introduction of cumbersome items

PERSONALIZE YOUR REELS ON REQUEST:

- Medical Paper 70 g/m² and 80 g/m² available
- Film color available also in Green and Transparent
- Add the FORMALDEHYDE indicator
- Flat Reels available also in length 100 m
- Print your Private Label on the items

FLAT REELS

Code	measure
RF0520	50mm x 200m
RF7520	75mm x 200m
RF1020	100mm x 200m
RF1220	120mm x 200m
RF1520	150mm x 200m
RF2020	200mm x 200m
RF2520	250mm x 200m
RF2720	270mm x 200m
RF3020	300mm x 200m
RF3220	320mm x 200m
RF3520	350mm x 200m
RF3820	380mm x 200m
RF4020	400mm x 200m
RF4220	420mm x 200m
RF5020	500mm x 200m
RF6020	600mm x 200m

GUSSETED REELS

Code	measure
RS7510	75mm x 25mm x 100m
RS1010	100mm x 50mm x 100m
RS1510	150mm x 50mm x 100m
RS2010	200mm x 50mm x 100m
RS2510	250mm x 65mm x 100m
RS3010	300mm x 65mm x 100m
RS3510	350mm x 65mm x 100m
RS4010	400mm x 65mm x 100m